



## MiniCollect® LH Lithium Heparin (Sep) Tubes



For In Vitro Diagnostic Use

### Instructions for use

#### Intended Use

To collect, transport, separate and process capillary blood for testing plasma in the clinical laboratory.

#### Product Description

**MiniCollect®** Tubes are plastic, non-evacuated, non-sterile low sample volume tubes with integrated collection devices. The closure is colour coded to identify the additives which are present in varying concentrations depending on the tube type and stated volumes.

The interior of the tube wall is coated with lithium heparin. The heparin anticoagulant activates antithrombin, thus blocking the coagulation cascade and producing a whole blood / plasma sample making it ideal for rapid analysis and analysis of blood from patients on anticoagulant therapy.

**MiniCollect®** Lithium Heparin Sep Tubes contain a gel on the bottom of the tube which moves upwards during centrifugation to form a stable barrier between the cells and the plasma.

The product is to be used by appropriately trained healthcare professionals in accordance with these instructions.

Tube Type	Matrix	Cap colour
<b>MiniCollect®</b> Lithium Heparin Tubes – not available in USA	Capillary blood	Green
<b>MiniCollect®</b> Lithium Heparin Sep Tubes (with gel) – not available in USA	Capillary blood	Mint green

#### Product versions

**MiniCollect®** Tubes with optional 13x75 mm carrier tubes (clear, amber)

**MiniCollect®** Complete Tubes pre-assembled with 13x75 mm carrier tubes

#### Storage before use

Store tubes at 4–25°C (40–77°F). Avoid exposure to direct sunlight. Exceeding the maximum recommended storage temperature may lead to impairment of the tube quality (*i.e. colouring, etc.*).

#### Sample Stability and Storage

Refer to the instrument assay instructions for use or literature for information on the correct sample material, correct storage and stability.

Heparin plasma should be separated from cells within 2 hours, either by collection and centrifugation with a gel tube or by transferring plasma into a secondary container if a gel tubes is not used.

#### Precautions/Cautions

- Insufficient or delayed mixing of tubes with additives may result in platelet clumping, clotting and/or incorrect test results.
- Do not use tubes if foreign matter is present.
- Handle all biological samples and blood collection devices according to the policies and procedures of your facility.
- Obtain appropriate medical attention in the case of any exposure to biological samples.
- Check all tubes to verify appropriate product and shelf-life before use. Do not use tubes after the expiration date.
- It is the laboratory's ultimate responsibility to verify that a change from one tube to another does not significantly affect analytical results obtained from patient samples.
- Single use only.
- To prevent contamination or leakage, make sure to avoid any blood in the recess around the scoop in the upper tube rim before closing the tubes. In case there is blood in the recess, it is recommended to clean it before closing the tube.

Only applicable for member states of the European Union: Should any serious incidents occur in relation to the product, these must be reported to the manufacturer and the competent authority in the member state, in which the user/patient is established.

#### Specimen Collection and Handling

##### Equipment required but not provided

- Labels for positive patient identification of samples
- For testing bilirubin, **MiniCollect®** tubes must be used with the amber carrier tube or other protective cover.
- Gloves and appropriate apparel for protection against exposure to blood borne pathogens
- Alcohol swab
- Dry gauze
- Warming device if appropriate, depending on the volume of blood needed and the test(s) to be performed
- Adhesive plaster or bandage
- Biohazard disposal container
- Safety lancet

##### Recommended Order of Draw (based on CLSI GP42-A6)

If multiple specimens are to be collected, including EDTA specimens, the EDTA specimen is drawn first to ensure adequate volume and accurate haematology test results. Specimens with other additives are collected next; serum specimens are collected last.

- 1 EDTA
- 2 Heparin / Heparin Sep
- 3 Glycolytic inhibitor tubes
- 4 Serum / Serum Sep

### Cap Removal

The ribbed area around the tube top indicates the cap opening location. Remove the cap by applying gentle pressure on the cap in an upward direction. The triangle indicator is positioned the opposite of the collection scoop.

### Specimen Collection

For the collection of blood by skin puncture, please refer to the instructions for use for the lancet device used.

Blood flow from the skin puncture site is increased by holding the puncture site in a downward position. The first drop of blood should be wiped away with a gauze pad if appropriate for the test method, as it is likely to contain excess tissue fluids. Apply gentle, intermittent pressure to the surrounding tissue (or proximal to the puncture site when using a finger). Strong, repetitive pressure (milking), must not be applied. It may cause hemolysis or tissue-fluid contamination of the specimen and adversely affect test results.

Drops of blood should be allowed to flow freely into the tube and down the walls of the **MiniCollect**<sup>®</sup> Tube. If a drop becomes lodged inside the scoop or to mix the contents as specimen is collected, gently tap the tube on a hard surface. Do not shake and avoid flicking the open tube. For correctly filled tubes, observe the fill mark. After collection, close the tube with the original cap, an audible click indicates correct closure. Invert until the blood completely mixes with the additive (approximately 8-10x). Vigorous shaking may cause foaming and hemolysis. The patient must be positively identified and the patient's blood sample must be properly labeled at the time of collection. The specimen should be labelled immediately following collection and mixing.





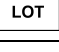

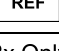
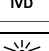


### Centrifugation

Ensure that tubes are properly seated in the centrifuge carrier. **MiniCollect**<sup>®</sup> tubes are recommended to be centrifuged at 3000g for a period of 10 minutes. Other centrifugation settings need to be validated by the user. Centrifugation should be done at a temperature of 15-24°C (59-75°F). It is not recommended to re-centrifuge separator tubes once the barrier has been formed. Barriers are more stable when tubes are spun in centrifuges with horizontal swing-out rotors rather than those with fixed angle heads.

### Disposal


- The general hygiene guidelines and legal regulations for the proper disposal of infectious material should be considered and followed.
- Always wear gloves during blood collection and disposal.
- Contaminated or filled blood collection tubes must be disposed of in suitable biohazard disposal containers, which can then be autoclaved and incinerated afterwards.

### Label Information

	Manufacturer		Temperature limit
	Use-by date		Do not re-use
	Batch code		Consult instructions for use
	Catalogue number		<i>In vitro</i> diagnostic medical device
	Prescription device (USA)		Keep away from sunlight

### Literature

CLSI. *Collection of Capillary Blood Specimens*. 7th ed. CLSI standard GP42. Clinical and Laboratory Standards Institute; 2020.

 Greiner Bio-One GmbH  
Bad Haller Str. 32  
4550 Kremsmünster, Austria

[www.gbo.com/preanalytics](http://www.gbo.com/preanalytics)  
[office@at.gbo.com](mailto:office@at.gbo.com)  
Phone +43 7583 6791

### Appendix – Log of changes since last revision

Location in the document	Change description
Precautions/Cautions	Addition of a precaution regarding blood in the recess around the scoop of the tube.
Precautions/Cautions	Addition of a note regarding reporting of serious incidences in the European Union.
Centrifugation	Clarification that other centrifugation settings must be validated by the user.
Label Information	Adaption of the <i>Rx Only</i> symbol to <i>Rx Only (USA)</i> .
Label Information	Addition of the <i>Keep away from sunlight</i> symbol.
Literature	Update of CLSI reference.