

Customer Information

Greiner Bio-One Multipurpose Beakers



The new Greiner Bio-One Multipurpose Beaker is manufactured from polypropylene. The beaker features an INTEGRITY SEAL, which ensures sterility and intactness of the product.

The beaker is leak-proof, shatterproof and resistant against chemical influences.

The advantages of this hygienic all-rounder:

- The INTEGRITY SEAL guarantees sterility - it is no longer necessary to individually pack each beaker
- Greiner Bio-One's conscious decision to reduce packaging materials is a contribution to environmental protection
- Save on disposal costs due to reduced amount of waste
- Timesaving, as there is no need to remove the beaker from a pack
- More information and better traceability - item number, lot number and expiry date are printed on the base of the beaker
- Writing space and exact graduation marks provide clarity



P. T. O.!


Customer Information

Greiner Bio-One Multipurpose Beakers

Handling Multipurpose Beakers with pre-assembled ring screw cap

(724402, 724461, 724473)

- 1) The ring screw cap is assembled on the beaker. To open twist off the cap, thereby perforating the INTEGRITY SEAL.
- 2) The ring then breaks off from the cap but remains on the beaker. Once opened the beaker is no longer sterile.

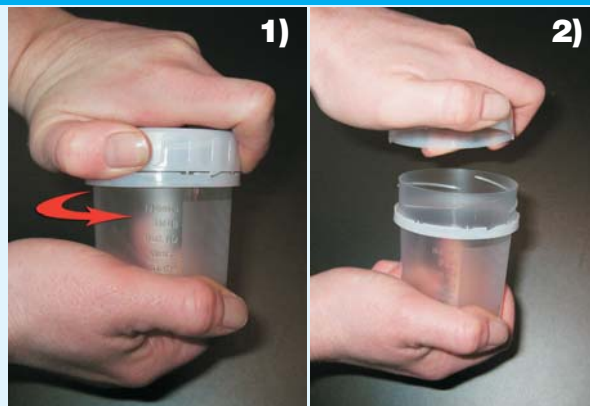

Handling Multipurpose Beakers without pre-assembled ring screw cap

(724401)

- 1) To apply the ring screw cap to the beaker, a light resistance must be overcome. The pressure is necessary to ensure that the integrity seal of the cap snaps onto the beaker.
- 2) When the beaker is opened, the rings break off from the cap but remain on the beaker.

Please note that the item numbers have changed.

The new beaker is available in the following item numbers:

Item n°.	Description	Packaging
724401	Multipurpose beaker, 100ml, without cap, graduated	1000
724170	Ring screw cap, white, for beaker (724401)	500
724173	Ring screw cap, red, for beaker (724401)	500
724402	Multipurpose beaker, 100ml, with ring screw cap, white, graduated, non-sterile	300
724461	Multipurpose beaker, 100ml, with ring screw cap, white, graduated, sterile	300
724473	Multipurpose beaker, 100ml, with ring screw cap, red, graduated, sterile	300