

FOR DIAGNOSTICS & IMMUNOLOGY

FREE OF
detectable
DNase

FREE OF
detectable
human DNA

FREE OF
detectable
RNase

 non-
cytotoxic

 non-
pyrogenic


96 WELL ELISA SINGLE-BREAK STRIP PLATES

Greiner Bio-One is offering C8 single-break strip plates with and without colour coding for diagnostic and immunological research applications. C-bottom single-break strips are supplied as twelve 8-well strips in a support frame with 96 spaces (12 x 8 matrix). The individual wells can be broken off separately ("single-break" option) and the number of tests performed can thus be precisely adjusted to the number of samples. The "C" describes a flat-bottom profile with rounded corners. The flat bottom enables precise optical measurements, whereas the rounded corners ensure that the individual wells can be pipetted without leaving a residue.

The quality of our immunological products is constantly checked in our quality assurance laboratory by means of ELISA. Clear MICROLON products are tested in a colorimetric immunoassay, black FLUOTRAC plates in a fluorescence immunoassay, white LUMITRAC plates in a luminescence immunoassay.

KEY FACTS

- / With medium and high-binding surfaces
- / Breakable C8 strips (clear) mounted in frame
- / With or without colour-coded well rims (blue / green / red / yellow)

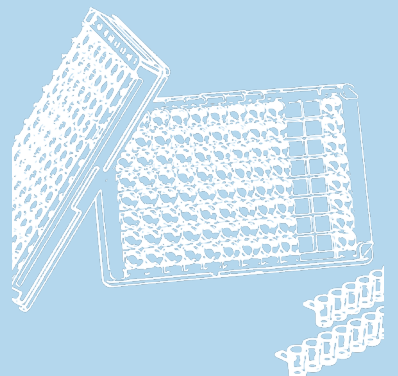

PRODUCT OVERVIEW

FOR FURTHER INFORMATION AND/OR SAMPLE ORDERING
PLEASE VISIT OUR WEBSITE OR CONTACT US.

Single-break Strip Plates

Well format: 96, Well profile: C-bottom, Bottom: solid, Raw material: PS, Plate design: strip plate 12x8, Working volume (well): 20 µl - 300 µl, Lid: no

Item no.	Binding characterisitc	Binding brand name	Colour stripe	Colour coding well rim	Qty. inner / outer
705070	med.-binding	MICROLON 200	○ clear		10 / 100
705070-25	med.-binding	MICROLON 200	○ clear		25 / 150
705071	high-binding	MICROLON 600	○ clear		10 / 100
705071-25	high-binding	MICROLON 600	○ clear		25 / 150
705063	med.-binding	MICROLON 200	○ clear	● red	10 / 100
705073	high-binding	MICROLON 600	○ clear	● red	10 / 100
705074	high-binding	MICROLON 600	○ clear	● blue	10 / 100
705065	med.-binding	MICROLON 200	○ clear	● green	10 / 100
705075	high-binding	MICROLON 600	○ clear	● green	10 / 100
705066	med.-binding	MICROLON 200	○ clear	● yellow	10 / 100
705076	high-binding	MICROLON 600	○ clear	● yellow	10 / 100