





SIMPLIFY YOUR CELL CULTURE WORK PROCESS

FREE OF detectable DNase	FREE OF detectable RNase	FREE OF detectable human DNA	FREE OF heavy metal
STERILE		 non-cytotoxic	 non-pyrogenic

CELLdisc

With closed filling caps (-CF1/-CF2)

Working in an GMP surrounding requires maximal security on sterility and hence a closed system for fluid transfer which excludes any opening of a cell culture disposable during cell cultivation. For such applications Greiner Bio-One has developed the new CELLdisc versions with closed filling caps. These consist of flexible USP class VI certified silicone tubes attached to the CELLdisc screw cap which can be safely connected to liquid reservoirs like media bags using the included MPC connector. The closed filling caps are preassembled with the respective CELLdisc and the entire disposable is sterilized as one unit, assuring maximal sterility, system integrity and highest customer convenience from the first moment.

For each CELLdisc two versions are available:

- / -CF1: closed filling cap with single tubing
- / -CF2: closed filling cap with double tubing with included dip-in tube

EXTERNAL FILTERS: FOR HIGHEST FILTRATION EFFICIENCY

CELLdisc closed filling cap variants are equipped with external filters of a defined pore size to maintain cell culture sterility and for ideal air flow. Thus, cells find the ideal growth conditions.



PRODUCT OVERVIEW

FOR FURTHER INFORMATION AND/OR SAMPLE ORDERING
PLEASE VISIT WWW.GBO.COM OR CONTACT US.

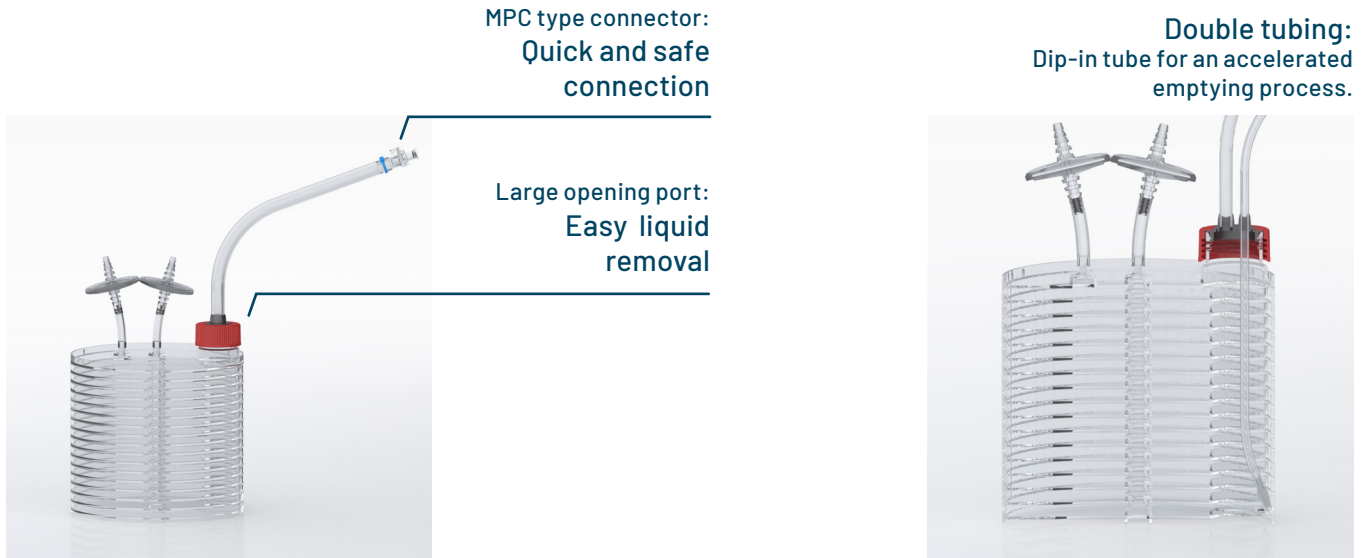


Figure 1: CELLdisc 16 layers with closed filling cap (-CF1) and CELLdisc 16 layers with closed filling cap (-CF2) and dip-in tube



Before first usage, the instructions of use must be read.

- / Download the instructions for use on our website www.gbo.com/downloads
- / Watch the handling video by scanning the QR-Code



CELLdisc with closed filling caps

Sterility SAL 10⁻⁶, triple packed

● Adv. TC on demand

Item No.	Layers	Surface	Description
6781xx-CF1	1-40	● TC	CELLdisc with closed filling cap, single tubing and external filter
6781xx-CF2	4-40	● TC	CELLdisc with closed filling cap, double tubing and external filter

CELLdisc with external filter

Sterility SAL 10⁻⁶, triple packed

● Adv. TC on demand

Item No.	Layers	Surface	Description
6781xx-EXF	1-40	● TC	CELLdisc with external filter

Intended Use: General laboratory products for cell culture to be used by qualified personnel in a laboratory environment.