

Your **Power** for Health



greiner bio-one



VACUETTE®

SAFETY Blood Collection Sets

Our Innovations
for Your Safety

www.gbo.com/preanalytics


SAFETY Blood Collection Sets

Specially for difficult vein conditions

The SAFETY Blood Collection Set has been specially developed for patients with difficult vein conditions , but is equally suitable for any regular blood taking situation, due to the gentle comfort it provides to a patient.

As the protection mechanism is activated whilst the needle is still in the vein, a high level of safety is ensured. Correct activation via retraction is indicated by an acoustic signal.

Protect yourself and your employees from the risk of needlestick injury. Every needlestick injury is one too many.

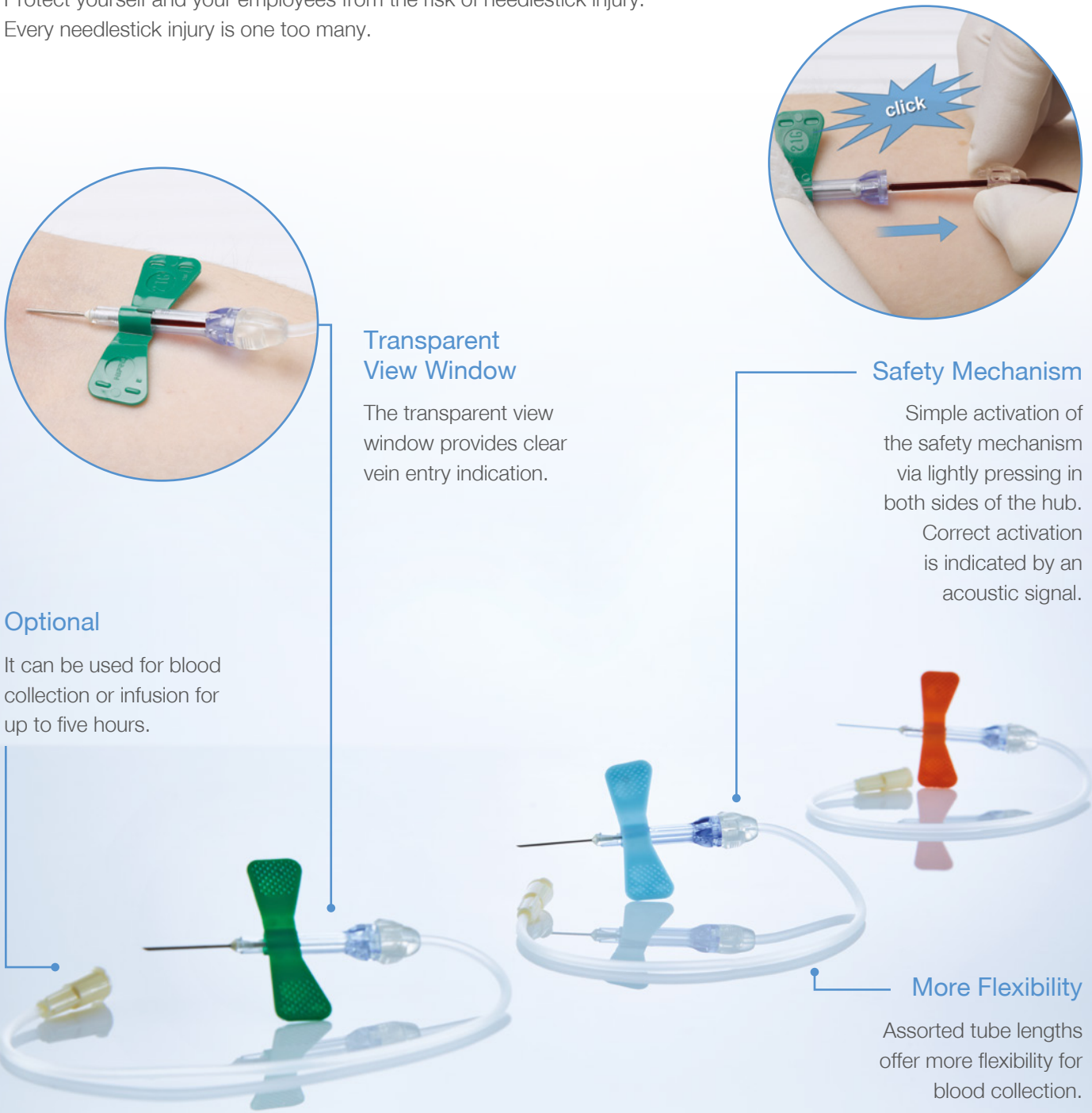

Transparent View Window

The transparent view window provides clear vein entry indication.

Safety Mechanism

Simple activation of the safety mechanism via lightly pressing in both sides of the hub. Correct activation is indicated by an acoustic signal.

More Flexibility

Assorted tube lengths offer more flexibility for blood collection.

SAFETY Blood Collection Sets Combination Products

As a practical combination product, the SAFETY Blood Collection Set is also available pre-assembled.

SAFETY Blood Collection Set + Luer Adapter

- Pre-attached set with Luer adapter in sterile blister pack.
- Handling is optimised.



SAFETY Blood Collection Set + Holder

- Pre-attached set with holder in sterile blister pack.
- No threading necessary.
- A practical, ready-to-use solution that saves time.

Item No. 450080
Item No. 450085-88
Item No. 450160-61

SAFETY Blood Collection Set + Blood Culture Holder

- Pre-attached set with Blood Culture Holder in sterile blister pack.
- Designed for use with most commonly used blood culture bottles and **VACUETTE®** blood collection tubes in one procedure.
- The set is sterile and pre-assembled for convenience, cost effectiveness and ensures your patient only has to undergo one venipuncture for all procedures.



Instructions for Use *



1 Remove protective cap




2 Venipuncture



3 Insert the tube



4 Remove the last tube

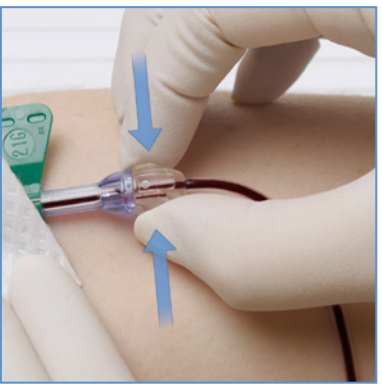

5 Activate the safety mechanism (by pressing down and pulling back the lighter coloured area)



6 Audible click indicates activation



7 Activated SAFETY Blood Collection Set



8 Disposal

* For the complete handling recommendation please refer to the Instructions for Use for each product.

SAFETY Blood Collection Sets

Item. No.	Description	Size		Packaging Inner Outer	
SAFETY Blood Collection/Infusion Sets (sterile)					
450090	SAFETY Blood Collection/Infusion Set	tubing length: 19cm	25G x 3/4"	50 pcs.	1000 pcs.
450091	SAFETY Blood Collection/Infusion Set	tubing length: 19cm	21G x 3/4"	50 pcs.	1000 pcs.
450191	SAFETY Blood Collection/Infusion Set	tubing length: 30cm	21G x 3/4"	50 pcs.	1000 pcs.
450092	SAFETY Blood Collection/Infusion Set	tubing length: 19cm	23G x 3/4"	50 pcs.	1000 pcs.
450192	SAFETY Blood Collection/Infusion Set	tubing length: 30cm	23G x 3/4"	50 pcs.	1000 pcs.
SAFETY Blood Collection Sets + Luer Adapters (sterile)					
450083	SAFETY Blood Collection Set + Luer Adapter	tubing length: 10cm	21G x 3/4"	50 pcs.	1000 pcs.
450081	SAFETY Blood Collection Set + Luer Adapter	tubing length: 19cm	21G x 3/4"	50 pcs.	1000 pcs.
450095	SAFETY Blood Collection Set + Luer Adapter	tubing length: 30cm	21G x 3/4"	50 pcs.	1000 pcs.
450084	SAFETY Blood Collection Set + Luer Adapter	tubing length: 10cm	23G x 3/4"	50 pcs.	1000 pcs.
450082	SAFETY Blood Collection Set + Luer Adapter	tubing length: 19cm	23G x 3/4"	50 pcs.	1000 pcs.
450096	SAFETY Blood Collection Set + Luer Adapter	tubing length: 30cm	23G x 3/4"	50 pcs.	1000 pcs.
450099	SAFETY Blood Collection Set + Luer Adapter	tubing length: 30cm	25G x 3/4"	50 pcs.	1000 pcs.
SAFETY Blood Collection Sets + Holders (sterile)					
450087	SAFETY Blood Collection Set + Holder	tubing length: 10cm	21G x 3/4"	24 pcs.	240 pcs.
450085	SAFETY Blood Collection Set + Holder	tubing length: 19cm	21G x 3/4"	24 pcs.	240 pcs.
450160	SAFETY Blood Collection Set + Holder	tubing length: 30cm	21G x 3/4"	24 pcs.	240 pcs.
450088	SAFETY Blood Collection Set + Holder	tubing length: 10cm	23G x 3/4"	24 pcs.	240 pcs.
450086	SAFETY Blood Collection Set + Holder	tubing length: 19cm	23G x 3/4"	24 pcs.	240 pcs.
450161	SAFETY Blood Collection Set + Holder	tubing length: 30cm	23G x 3/4"	24 pcs.	240 pcs.
450080	SAFETY Blood Collection Set + Holder	tubing length: 30cm	25G x 3/4"	24 pcs.	240 pcs.
SAFETY Blood Collection Sets + Blood Culture Holders (sterile)					
450182	SAFETY Blood Collection Set + Blood Culture Holder	tubing length: 19cm	21G x 3/4"	24 pcs.	240 pcs.
450183	SAFETY Blood Collection Set + Blood Culture Holder	tubing length: 19cm	23G x 3/4"	24 pcs.	240 pcs.
450184	SAFETY Blood Collection Set + Blood Culture Holder	tubing length: 30cm	21G x 3/4"	24 pcs.	240 pcs.
450185	SAFETY Blood Collection Set + Blood Culture Holder	tubing length: 30cm	23G x 3/4"	24 pcs.	240 pcs.

These products are not made with natural rubber latex or DEHP.



Austria (Main Office)

Greiner Bio-One GmbH
Phone +43 7583 6791-0
E-Mail office@at.gbo.com

Hungary

Greiner Bio-One Hungary Kft.
Phone +36 96 213 088
E-Mail office@hu.gbo.com

Brazil

Greiner Bio-One Brasil
Phone +55 19 3468-9600
E-Mail office@br.gbo.com

India

Greiner Bio-One INDIA Pvt., Ltd.
Phone +91 120 375 9291
E-Mail info@gboindia.com

Netherlands

Greiner Bio-One B.V.
Phone +31 172 4209 00
E-Mail info@nl.gbo.com

Thailand

Greiner Bio-One Thailand Ltd.
Phone +66 38 4656 33
E-Mail office@th.gbo.com

China

Greiner Bio-One Suns Co., Ltd.
Phone +86 10 83 55 19 91
E-Mail office@cn.gbo.com

Indonesia

Greiner Bio One Indonesia KPPA
Phone +62 812 1391 3900
E-Mail office@id.gbo.com

Portugal

VACUETTE Portugal S.A.
Phone +351 252 647 721
E-Mail info@vacuette.pt

Turkey

Greiner Bio-One Turkey
Phone +90 216 576 6004
E-Mail office@tr.gbo.com

France

Greiner Bio-One SAS
Phone +33 1 69 86 25 25
E-Mail accueil@gbo.com

Italy

Greiner Bio-One Italia S.r.l.
Phone +39 02 9438 3340
E-Mail office@it.gbo.com

Spain

Greiner Bio One España, S.A.U.
Phone +34 91 652 77 07
E-Mail info@es.gbo.com

United Kingdom

Greiner Bio-One Ltd.
Phone +44 1453 8252 55
E-Mail info@uk.gbo.com

Germany

Greiner Bio-One GmbH/Prenalytics
Phone +49 7022 948-0
E-Mail office@de.gbo.com

Japan

Greiner Bio-One Co. Ltd.
Phone +81 3 3505 8875
E-Mail info@jp.gbo.com

Switzerland

Greiner Bio-One VACUETTE
Schweiz GmbH
Phone +41 7 12 28 55 22
E-Mail office@ch.gbo.com

USA

Greiner Bio-One North America Inc.
Phone +1 704 261-7800
E-Mail info@us.gbo.com