



VEIN VISUALIZATION IN A NEW LIGHT


CHKISTIE+

VEINVIEWER® Flex

Increasingly more patients with difficult vein conditions urgently need a solution. VeinViewer® Flex is the solution to meet these needs.

VeinViewer® Flex combines all the recent advances in the visualization of veins into one portable device. VeinViewer® Flex helps to find veins more quickly, especially for problem patients, making puncture much easier. The Flex device delivers a revolutionary, real-time HD image that visualizes the vein trajectories before, during, and after puncture. The device can be used to display peripheral veins up to 10 mm below the skin or blood patterns up to 15 mm below the skin surface.

- / Fast battery swap, long run time.
- / Reliable in daily routines.
- / Sturdy S-mount is quickly assembled.
- / Ergonomic design for ease of use.
- / Featuring the integrated graphic catheter size reference card, CathCompass™

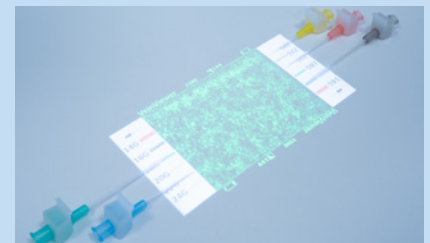

ASSESS™ IMAGING SUITE.

See veins like never before, in real-time up to 10mm below the skin. A wide range of different modes allows for optimal vein visualization depending on the individual patient and the lighting conditions. The integrated graphic reference card CathCompass helps choose the right sized IV cannula.

FOR MAXIMUM CONTROL:

- / Fine details for small veins
- / Image acquisition is HIPAA compliant
- / Inverse colour mode
- / Size adjustment in three different window sizes
- / MaxBright for changing lighting conditions
- / CathCompass™ for improved decision making



PRODUCT OVERVIEW

FOR FURTHER INFORMATION AND/OR SAMPLE ORDERING
PLEASE VISIT OUR WEBSITE OR CONTACT US.

VeinViewer® Flex

Item No.	Description	Inner / Outer [Qty.]
890001	Provided with 2 pc. lithium batteries, a charging station, 2 flexible mounting arms and 2 clamps	-/1

VeinViewer® Accessories

Item No.	Description	Inner / Outer [Qty.]
890011	Charging station for rechargeable battery for VeinViewer® Flex	-/1
890011A	Adapter for charging station for rechargeable battery for VeinViewer® Flex	-/1
890010	Spare rechargeable battery for VeinViewer® Flex and Vision 2	-/1

RECOMMENDED HANDLING

Please refer to the user manual for detailed instructions



1

With two fast swap rechargeable batteries, there is no risk of downtime. Can also be powered by the mains if necessary.




2

The Flex can be easily positioned and operated with just one hand. If you need both hands free for the puncture, the S-mount quickly gives a hands-free approach.



3

The Flex is positioned parallel to the puncture site. The correct focus position can be identified by the fact that the text in the edge of the image is clearly legible (approx. 30 cm distance).

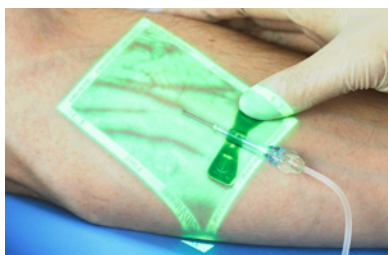

4

Standard mode also includes the CathCompass graphic, this can be deactivated by pressing F2.