

URINE COLLECTION SYSTEMS



VACUETTE®

Urine Collection Products

Urine Tubes

Z Urine No Additive Tubes are used as collection and/or transport tubes in chemical urine analysis, urine test strips and urine microscopy. The specimen should be refrigerated (2°C - 8°C) if it cannot be tested within 2 hours of collection.

The VACUETTE® Urine CCM Tubes contain a powdered stabiliser for bacterial and yeast culture. Urine samples collected in this tube can be stored at room temperature (20°C - 25°C) for up to 48h. The CCM tube can thus win time for carrying out the analysis.

Urine Accessories

The Urine Transfer Device allows easy and hygienic collection of sample material from a container into the tube. The product is available in two lengths.

Urine cups from Greiner Bio-One come with a transfer device already integrated in the cap. This feature gives you the advantage of being able to perform the sample transfer cleanly.

/ See next page for product and ordering details.

OPTIONS FOR SAFE SAMPLE TRANSFER



Urine Transfer Device for easy and hygienic sample transfer



Urine Cup with Integrated Transfer Device



PRODUCT OVERVIEW

FOR FURTHER INFORMATION AND/OR SAMPLE ORDERING
PLEASE VISIT OUR WEBSITE OR CONTACT US:

PHONE: **01453 825255**


EMAIL: **SALES.UK@GBO.COM**

VACUETTE® Z Urine No Additive Tubes



Item No.	NPC Code	Nominal volume	Cap colour	Cap type	Ring colour	Tube size	Additive	Pack
455007	KFK707	10 ml	yellow	pull	yellow	16 x 100	none	50
456007	KFK708	6 ml	yellow	premium (twist)	yellow	13 x 100	none	50

VACUETTE® Urine CCM Tubes

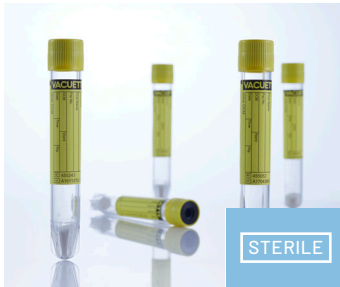


Item No.	NPC Code	Nominal volume	Cap colour	Cap type	Ring colour	Tube size	Additive	Pack
456284	KFK990	6 ml	yellow	pull	black	13 x 100	CCM (stabiliser)	50

Handling Video



VACUETTE® Urine Transfer Device

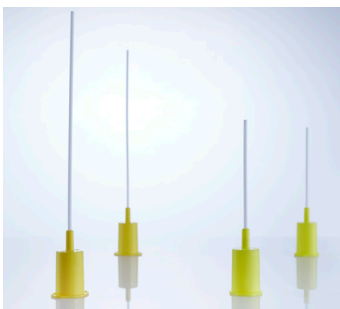

Item No.	NPC Code	Description	colour	Pack
450251	KFK435	VACUETTE® Urine transfer device, short version	yellow	50
450252	KFK569	VACUETTE® Urine transfer device, long version	yellow	50

VACUETTE® Urine Cup



Item No.	NPC Code	Description	Cap colour	Cap type	Pack
724321	-	Urine Cup with intergrated Transfer Device	yellow	screw cap	300

Handling Video

