


Your **Power** for Health



greiner bio-one

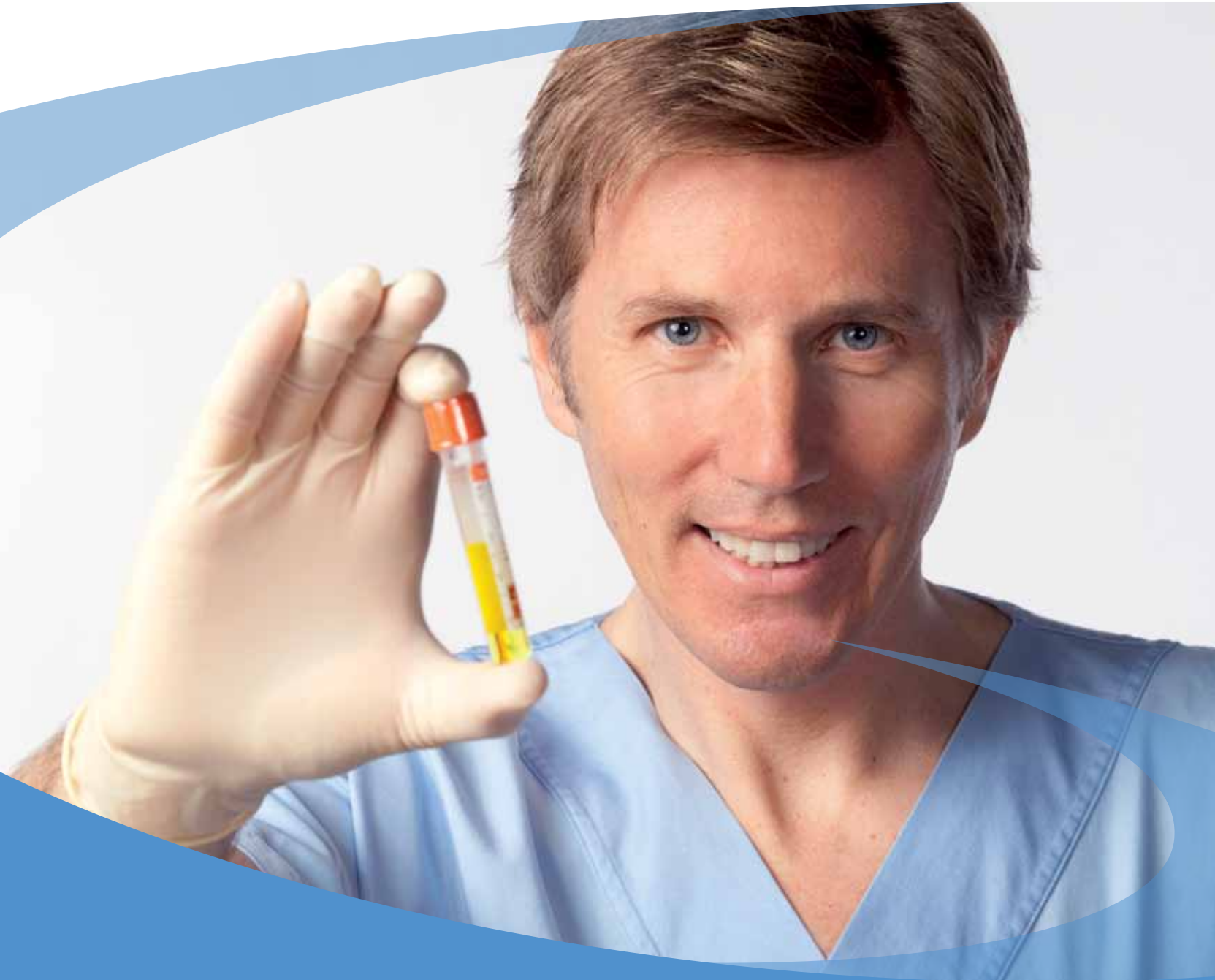

**SCS**

Saliva  
Collection System

## Saliva Collection System in saliva veritas

[www.gbo.com/preanalytics](http://www.gbo.com/preanalytics)

# GBO Saliva Collection System - SCS in saliva veritas



Greiner Bio-One has revolutionised saliva analytics with the new, unique collection and quantification system. This is the first system in the world on a liquid collection basis, and can ensure a representative collection from the entire oral cavity. The system fulfils the requirements of the In Vitro Diagnostics and Medical Device guidelines and provides application opportunities for developing new saliva diagnostics.

With this unique product, the photometric determination of the saliva quantity allows for the definite determination of the analytes in the initial volume.



## ... simple



The simplicity of the system means that the collection can be carried out by both patient and medical staff.

## ... safe



Protective film on the transfer tubes indicate if there has been unauthorised intervention. This ensures an easy visual indicator of any tampering to the sample or a substitution.

## ... fast



As a sufficiently sized sample can be collected at any time, delays are avoided.

# The Greiner Bio-One Saliva Collection System makes preanalytics even simpler

By rinsing with the saliva collection solution, saliva is collected easily and without stress. The dye contained (tartrazine\*) serves as internal standard and ensures reproducible saliva analysis. Sample collection with the collection beaker is then a comfortable procedure. The hygienic sample transfer is ensured due to the vacuum in the saliva transfer tube.

## The Components



## Product Overview

Item n°	Description	Packaging	
		Inner	Outer
881001	GBO Saliva Collection System, single-pack (for layperson)	1 pc.	24 pcs.
881002	GBO Saliva Collection System, 50 pcs. bulk Bulk-packed in a carton. All components are packed in one carton, sufficient for 50 applications.	1 carton	50 pcs.
881112	(1) GBO Saliva Extraction Solution 4ml	24 racks	1200 pcs.
881210	(2) GBO Saliva Collection Beaker, with integrated transfer unit	1 carton	200 pcs.
881113	(3) GBO Saliva Transfer Tube 3.5ml	24 racks	1200 pcs.

\*Tartrazine is permitted by the EU/FDA (directive 94/36/EC) for use as food colourant and can be found in various concentrations in cheese, spices, ice cream etc.



For further information, please visit our website [www.gbo.com/preanalytics](http://www.gbo.com/preanalytics) or contact us:

**Austria (Main Office)**  
Greiner Bio-One GmbH  
Phone +43 7583 6791-0  
Fax +43 7583 6318  
E-Mail [office@at.gbo.com](mailto:office@at.gbo.com)

**Brazil**  
Greiner Bio-One Brasil  
Phone +55 19 3468-9600  
Fax +55 19 3468-9621  
E-Mail [office@br.gbo.com](mailto:office@br.gbo.com)

**China**  
Greiner Bio-One Suns Co., Ltd.  
Phone +86 10 83 55 19 91  
Fax +86 10 63 56 69 00  
E-Mail [office@cn.gbo.com](mailto:office@cn.gbo.com)

**France**  
Greiner Bio-One SAS  
Phone +33 1 69 86 25 25  
Fax +33 1 69 86 25 35  
E-Mail [office@fr.gbo.com](mailto:office@fr.gbo.com)

**Germany**  
Greiner Bio-One GmbH/Preanalytics  
Tel +49 7022 948-0  
Fax +49 7022 948-514  
E-Mail [office@de.gbo.com](mailto:office@de.gbo.com)

**Hungary**  
Greiner Bio-One Hungary Kft.  
Phone +36 96 213 088  
Fax +36 96 213 198  
E-Mail [office@hu.gbo.com](mailto:office@hu.gbo.com)

**India**  
Greiner Bio-One INDIA Pvt., Ltd.  
Phone +91 120 456 8787  
Fax +91 120 456 8788  
E-Mail [info@gboindia.com](mailto:info@gboindia.com)

**Netherlands**  
Greiner Bio-One B.V.  
Phone +31 172 4209 00  
Fax +31 172 4438 01  
E-Mail [info@nl.gbo.com](mailto:info@nl.gbo.com)

**Spain**  
VACUETTE Espana S.A.  
Phone +34 91 652 77 07  
Fax +34 91 652 33 35  
E-Mail [info@vacuette.es](mailto:info@vacuette.es)

**Switzerland**  
Greiner Bio-One VACUETTE  
Schweiz GmbH  
Phone +41 7 12 28 55 22  
Fax +41 7 12 28 55 21  
E-Mail [office@ch.gbo.com](mailto:office@ch.gbo.com)

**Thailand**  
Greiner Bio-One Thailand Ltd  
Phone +66 38 4656 33  
Fax +66 38 4656 36  
E-Mail [office@th.gbo.com](mailto:office@th.gbo.com)

**United Kingdom**  
Greiner Bio-One Ltd.  
Phone +44 1453 8252 55  
Fax +44 1453 8262 66  
E-Mail [info@uk.gbo.com](mailto:info@uk.gbo.com)

**USA**  
Greiner Bio-One North America Inc.  
Phone +1 704 261-7800  
Fax +1 704 261-7899  
E-Mail [info@us.gbo.com](mailto:info@us.gbo.com)