

96 Well ELISA Single-Break Strip Plates

For Diagnostics and Immunological Research

Greiner Bio-One is offering C8 Single-break Strip plates with and without colour coding for diagnostic and immunological research applications.


C-bottom single-break strips are supplied as twelve 8-well strips in a support frame with 96 spaces (12 x 8 matrix). The individual wells can be broken off separately ("single-break" option) and the number of tests performed can thus be precisely adjusted to the number of samples.

The "C" describes a flat-bottom profile with rounded corners. The flat bottom enables precise optical measurements, whereas the rounded corners ensure that the individual wells can be pipetted without leaving a residue.

The quality of our immunological products is constantly checked in our quality assurance laboratory by means of ELISA. Clear MICROLON® products are tested in a colorimetric immunoassay, black FLUOTRAC™ plates in a fluorescence immunoassay, white LUMITRAC™ plates in a luminescence immunoassay.

Key Facts

- Strips mounted in frame
- C8 clear strips
- Optional red, blue, green or yellow colour coding on the well rim
- Medium and high binding surfaces
- L: 127.76 mm / W: 85.48 mm / H: 14.5 mm

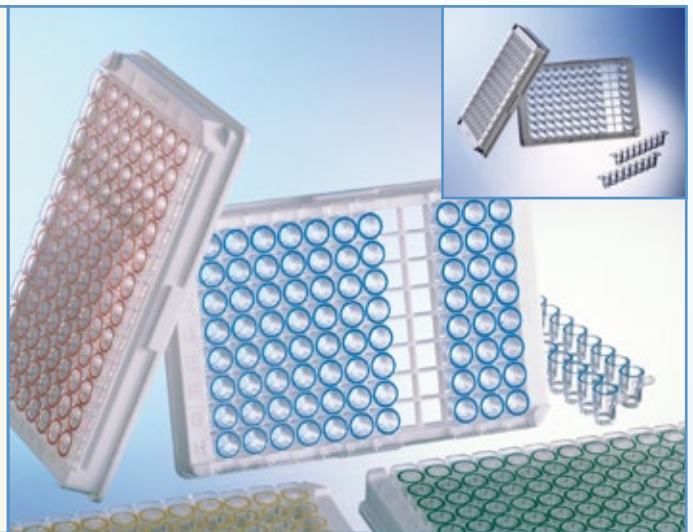

Please turn page for ordering information!

Ordering Information			
Cat. No.	Product Description	Quantity per Bag	Quantity per Case
705 070	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 200, med. binding, clear	10	100
705 071	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 600, high binding, clear	10	100
705 063	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 200, med. binding, clear, red colour coding	10	100
705 073	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 600, high binding, clear, red colour coding	10	100
705 074	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 600, high binding, clear, blue colour coding	10	100
705 065	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 200, med. binding, clear, green colour coding	10	100
705 075	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 600, high binding, clear, green colour coding	10	100
705 066	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 200, med. binding, clear, yellow colour coding	10	100
705 076	Single-break strip plate, 12 x C8 strips, PS, C-bottom, MICROLON® 600, high binding, clear, yellow colour coding	10	100